


GentloTM

RADIOFREQUENCY and PLASMA



Modular device

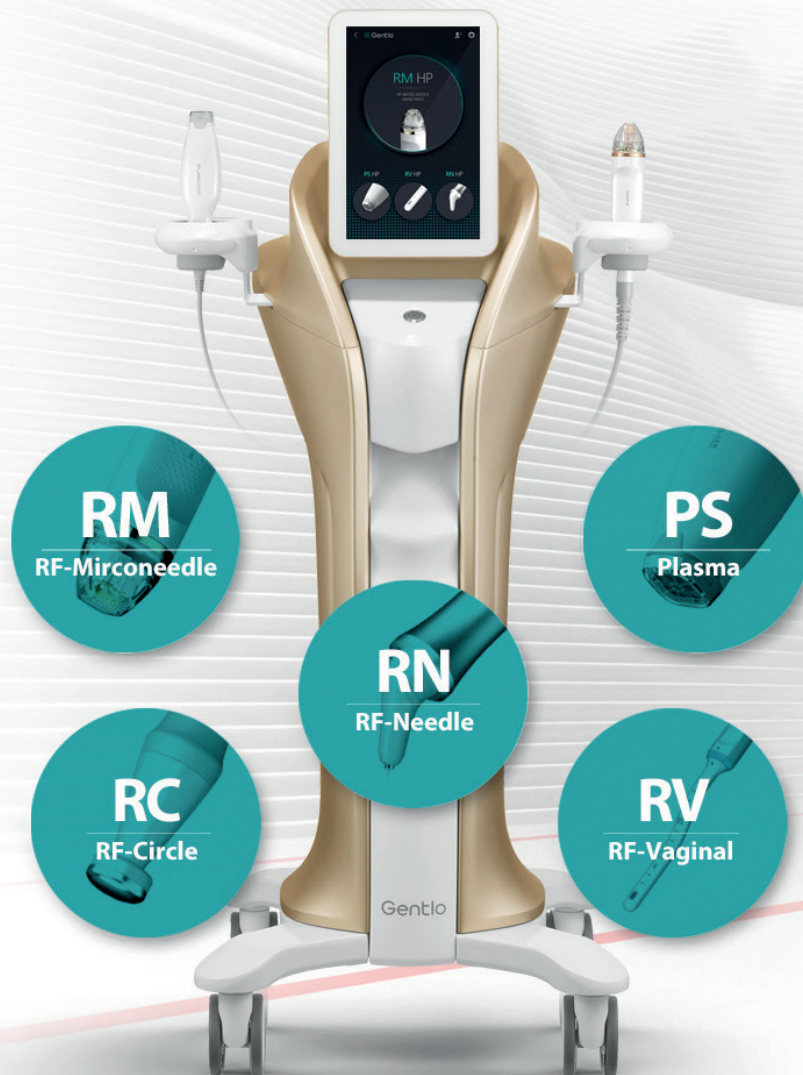


Good Design Selection
2019



Capactuel
MEDICAL

Preventive and Restorative Aesthetic

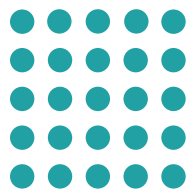

GentloTM

An Innovative combinaison of RF and PLASMA

GENTLO is the first device combining
RF and PLASMA technology.
Lifting effect, skin tightening,
overall skin revitalization,
from epidermis to dermis, an optimal
result thanks to a synergy of treatments.



RM | RF-Microneedle

Fractional RF microneedles



RN | RF-Needle

Monopolar RF Emitting Needle



PS | PLASMA

Plasma treatment head



RC | RF-Circle

Monopolar RF Flat Emitting Head



RV | RF-Vaginal

Monopolar RF Vaginal Probe




RC | RF-Circle

Monopolar RF Emitting Flat Head

Bulk Heating before and after treatment

With the small RF ceramic, the transmission of RF energy
on the skin is possible where it is difficult to use a needle.
To be used before and after acne treatment, overall skin rejuvenation.
Can be applied as part of a vulvar treatment.

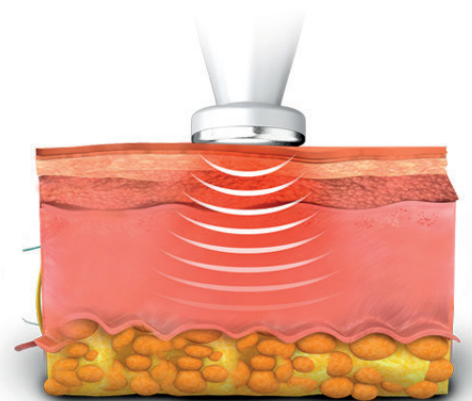


RV | RF-Vaginal

Monopolar RF Vaginal probe

Special vaginal treatment

Practical radiofrequency treatment for women.
Non-invasive, non-surgical, without anesthesia.
Vaginal rejuvenation and tightening, colpoxerosis,
it is possible to obtain an optimal treatment
thanks to an ergonomic curved handpiece.


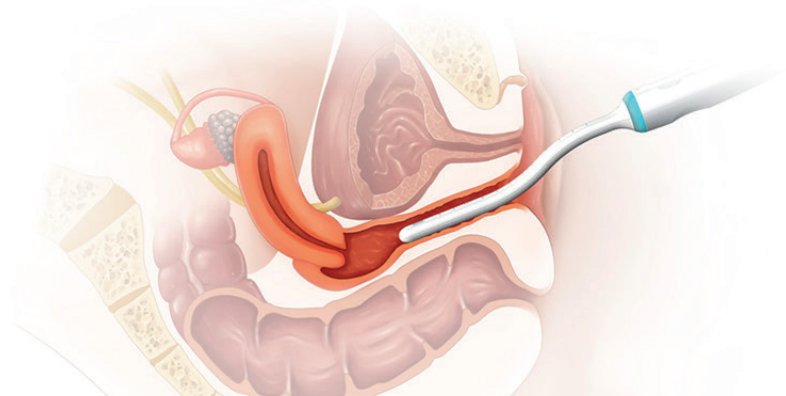




PS | PLASMA

Head for plasma treatment. Bleaching, sterilization, tightening, relaxation, absorption of cosmetic ingredients increased by 466%. Treatment results consolidated by clinical trial results.

Technology

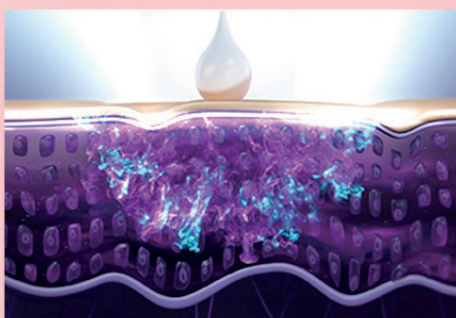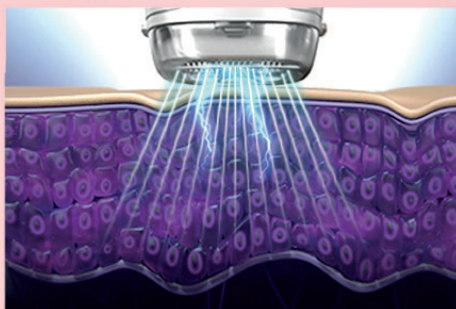

Generating high power and stable plasma.

Delivering plasma ion to the plasma tip.

Plasma affects the membrane potential for deep absorption of molecules solution.

Plasma sterilizes and stimulates regeneration of skin cells.



PLASONIC Skin Booster

The ampoule contains amino acids, vitamins and 4 different tripeptides that include a high concentration of growth factors. This is a double skin booster that improves the effect of PLASONIC treatments in skin rejuvenation and bleaching.




Moisturizing effect : 100%

Pores tightening : 82%

Softening effect : 91%

Wrinkles reduction : 91%

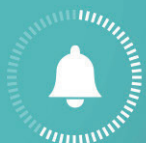

Alarm indicating the shots ending



Voice guidance

안녕하세요!
hello!
您好!
¡Hola!
أهلاً

5 languages available

AUTO

Selection mode



Memory function




RN | RF-Needle

Monopolar RF Emitting Needle


Gentlo™

Selective destruction of sebaceous glands

RF needle destroys selectively sebaceous glands.
Does not damage the epidermis, treats acne, prevents recurrence.

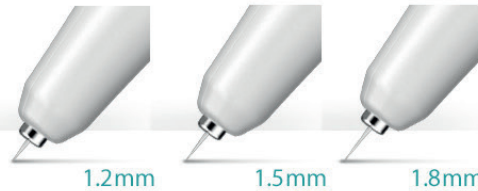

1.2mm, 1.5mm, 1.8mm (Isolated type) Thermic effect : 0.9mm, 1.2mm, 1.5mm


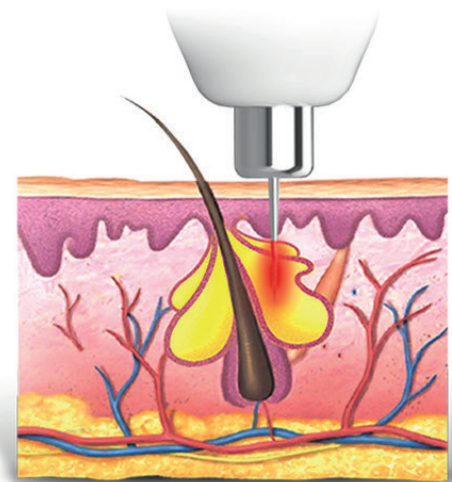


RM | RF-Microneedle

Microneedles for Fractional RF. Single and non-insulated needle for 2000 shots treatments.
Resurfacing, acne scar, hypertrophy, stretch mark, tightening.



GENTLO RF Special Micro Needles

Adjustable needle depth from 0.5 mm to 3.5 mm for insulated or insulated needle effect.

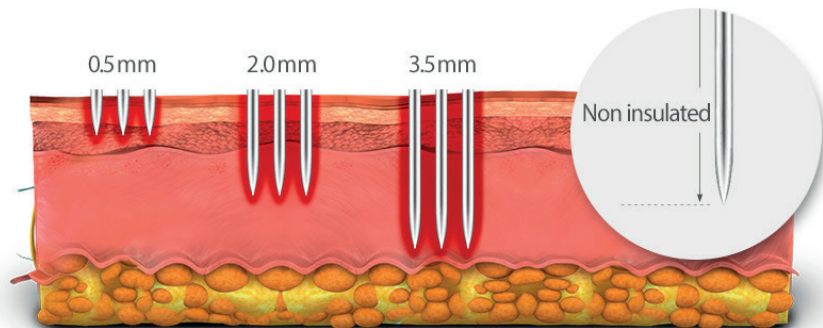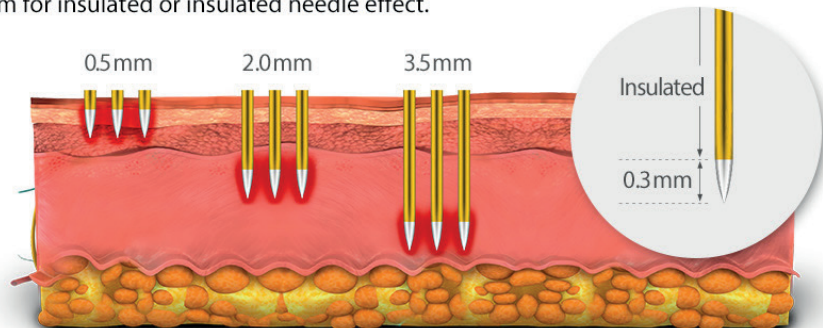

Wide view angle


25% Reduced tip surface size minimizes overlap during application that provides a wide viewing angle for a more focused and accurate procedure.

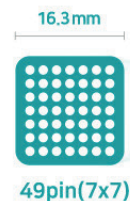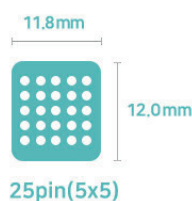

Lesions treatment

RF treatment

Treats scars, stretch marks, hyperhidrosis, pores, fine wrinkles, Allows overall skin rejuvenation, lifting, revitalizing, resurfacing, skin tightening.

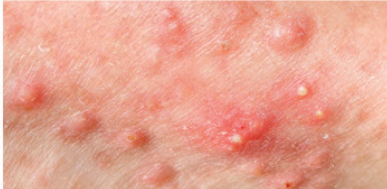

Plasma treatment

Improves texture and skin density and pore quality, Reduces wrinkles and dryness, Whitening, toning and moisturizing. Active ingredients administration.




RF for woman

Vaginal tightening, vulvar bleaching, vaginal rejuvenation.



1 Smart Integrated System (SIS)

Integrated intelligent system. Self-monitoring and online diagnostics for real-time troubleshooting.

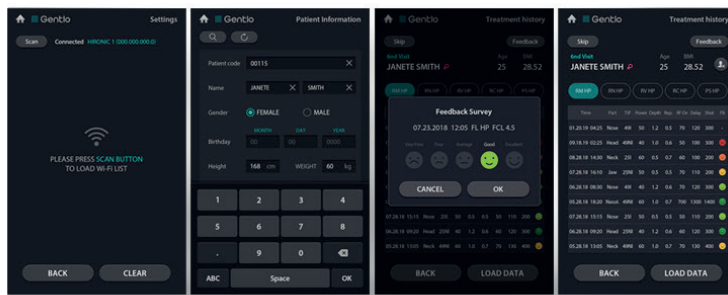

Remote management with Cloud System

- Notice of Replacement of Consumables
- Saved Protocols
- Self-diagnosis and inspection
- Prevention of malfunctions and support functions
- Checking the shots and the use of the camera
- Remote software/firmware update available


Search "HIRONIC" on the App store

2 Big Data System

Accumulation of processing data that allows GENTLO users to improve and secure protocols.

3 Graphical User Interface Simple & Intuitive

GENTLO Color User Interface

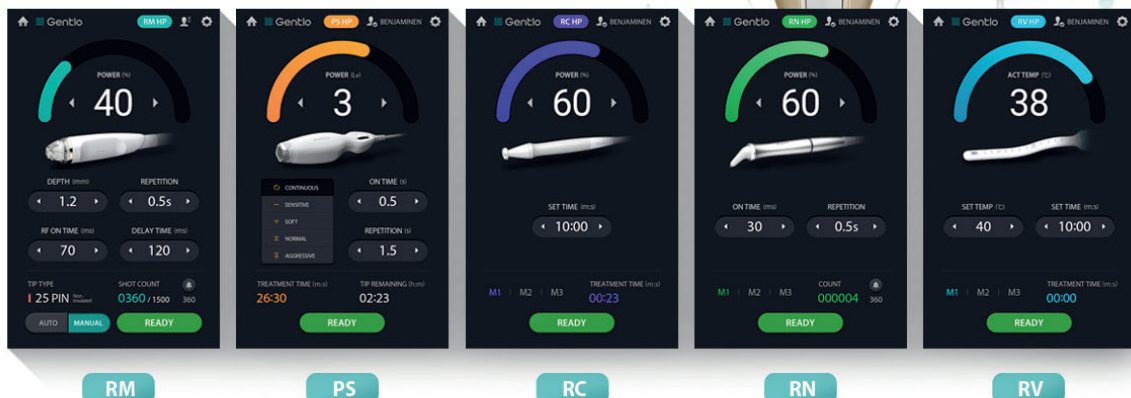

12.1" Display LCD

Maximum screen tilt of 75°



Sensor Technology

Handpiece recognition automatic system.



Ventilation system

Orifice for fast and powerful air circulation.



Handpiece	Radiofrequency	RF-Micro Needle	25pin	Insulated
				Non insulated
		49pin	Insulated	
			Non insulated	
		RF-Circle	30V~360V	
	RF-Needle 1pin	10V~168V		
	RF-Vaginal	10V~168V		
	PLASMA	Level 1~5		
Shot count		2000 shots		
Penetration depth		0.5~3.5mm (0.1 step)		
Repetition		Single / 0.2 / 0.5 / 0.8 / /1.0 / 2.0 (sec)		
RF - Intensity		Level 1~10		
Treating size		25PIN (8x8mm) / 49PIN (2x2mm)		
Frequency		2MHz		
RF emission time		50~950ms		
Convenience system		12.1" Wide LCD		
		5 Languages available (EN, CN, KR, ES, AR)		
		Voice guidance		
		Auto selection mode		
		Memory function		
		Air passing hole		
General	Max output	25W (25pin) / 50W (49pin)		
	Size	654(W) x 768(L) x 1247(H)		
	Weight	43kg		

Capactuel
MEDICAL

Preventive and Restorative Aesthetic