

MICCOL-S

Selective **Cryolipolysis** device
for natural body fat reduction.



- > 2 areas treated simultaneously.
- > Multi applicators.
- > Massage function.
- > Temperature control sensors.
- > 1 hour, it's done.

No surgery, No unavailability !
Painless, effective and safe.

www.capactuel.com

Capactuel
MEDICAL

Preventive and Restorative Aesthetic

MICCOOL-S

Cryolipolysis System

> WHAT IS CRYO FAT RÉDUCTION ?





MICCOOL-S is for people who wish to sustainably eliminate unsightly fat deposits that are resistant to diets and sport. Used Cryolipolysis technique is not a weight loss solution for obese people but a non-surgical alternative to liposuction.

MICCOOL-S requires no dietary supplements, drugs, or dietary plan to be effective. Once the results observed after treatment, a healthy lifestyle is self-evident to the patient.

> STRENGTHS OF MICCOOL-S CRYOLIPOLYSIS DEVICE.

- > Possibility of treating two areas simultaneously with two applicators. (Time saving, increased profitability).
- > 8 easily interchangeable nozzles tailored to specific areas to be treated.
- > Intuitive and ergonomic software to adjust : vacuum pressure, treatment time, temperature.
- > Each applicator can be set separately in order to optimize the results according to the targeted area.
- > Applicators have 8 temperature control sensors for greater efficiency and safety.
- > Vacuum massage function (better distribution of cells frozen by the cold, promotes the elimination of fat).

> TREATMENT PRINCIPLE.

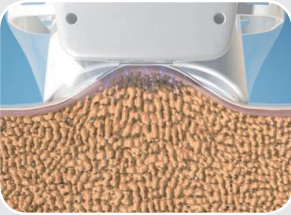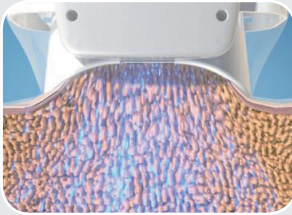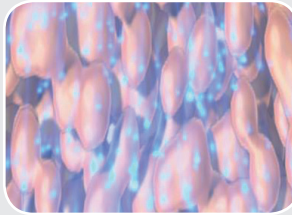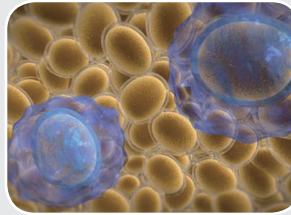			
Vacuum with controlled aspiration	Cold energy transmission	Fat cells cristallization*	Apoptosis or programmed cellular death

* Massage fonction at the end of treatment

> AN EFFICIENT AND NATURAL PROCESS.


MICCOOL-S causes apoptosis process which gradually reduces the layer of fat, although an immediate transformation is observed after the session. The removal of fat cells in the body is naturally carried out for 1 to 3 months beyond which the most remarkable results are usually observed.

> PRACTICAL, RELIABLE AND SAFE HAND PIECE.

MICOOL-S has 2 applicators with interchangeable tips for treating with greater ease, efficiency and safety different areas of the body.



Arms, belly, hip... In one hour, it is possible to treat 2 areas simultaneously !

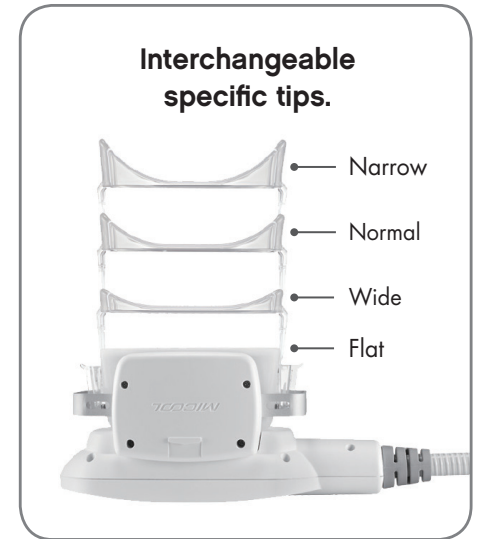


> A SAFE PROCEDURE.

Without surgery, MICOOL-S is non-invasive procedure without pain, or downtime. The cold energy is delivered to the subcutaneous fat layer by applicators without causing damage to surrounding tissues.

During the cryolipolysis process on the treated areas, 8 temperature control sensors, a Smart Gel Pad, a Smart Liner and a massage function ensure safety.


> APPLICATION ZONES.

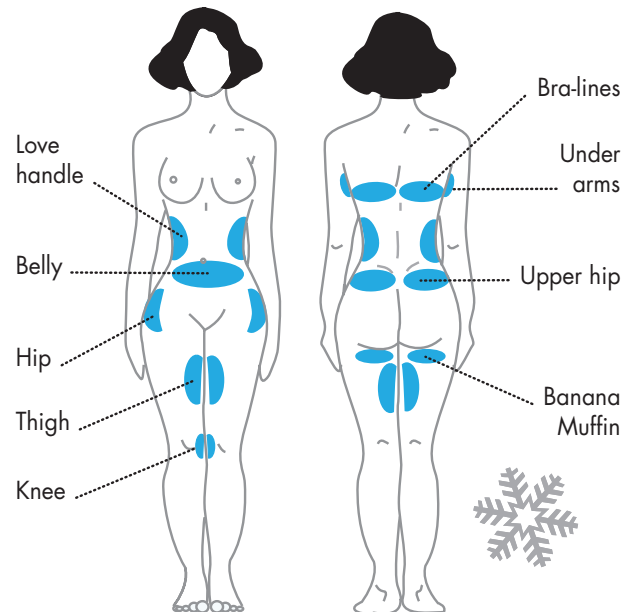

MICCOOL-S

Cryolipolysis System

> TREATMENT RESULT.

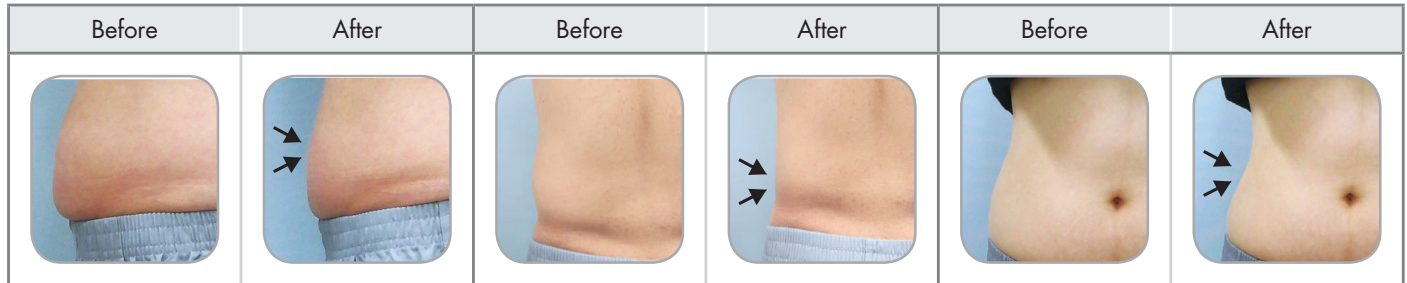


> SMART GEL PAD AND SMART LINER COATING.

During the energy transfer, the Smart Gel Pad protects the skin surface. It allows to use the correct temperature to get desired results even on the most sensitive area. Smart Liner coating contributes to personal hygiene, protects the machine against the absorption of foreign bodies, improves comfort suction and prolongs the life of the applicators..

> TECHNICAL FEATURES.

Energy type	Vacuum + Cooling + Massage
Treatment time	Smart Applicator : 30, 45, 60 mins
Treatment area size	Approximately 10 x 20 cm
Display	10.4 Inch LCD Touch Screen
Vacuum	- 310 mb
Electrical requirements	AC 100-220V, 50/60Hz
Dimensions (W x D x H)	340 x 470 x 1160 mm
Weight	63 kg

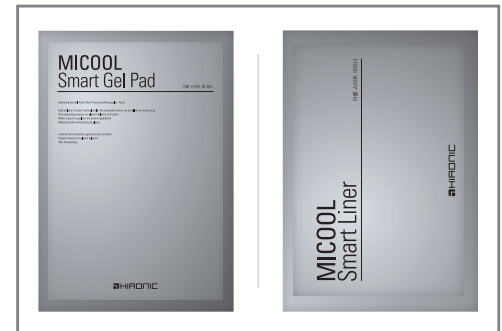

Smart Gel Pad and Smart Liner

Reveal your potential ! The secret of your success...

33 +(0)1 85 34 31 37 - info@capactuel.com - www.capactuel.com

Capactuel
MEDICAL

Preventive and Restorative Aesthetic